



HERVIK[®]
Healthcare Pvt Ltd



Broucher:
M/s. Hervik Healthcare Private Limited.

e-mail: hervikhealthcare@gmail.com
website: www.hervikhealthcare.com
Customer care: +91 9985 979 770.


Hypodermic Auto-disable (AD)/ Safety Disposable syringes



(Ribbon and Blister Packed “Higgs & Hervik” brands Hypodermic Auto-disable / safety Disposable Syringes with needle and without needle).



Auto-disable (AD) syringes automatically become disabled after one use. The "Higgs" syringe is a single-use, disposable, auto-disable syringe with a safety plunger that breaks off after a single use. The "Higgs" syringes can be purchased with a detachable needle. Auto disable syringes are meets ISO 7886-3 & ISO7886-4 Standards for A.D. Syringes.

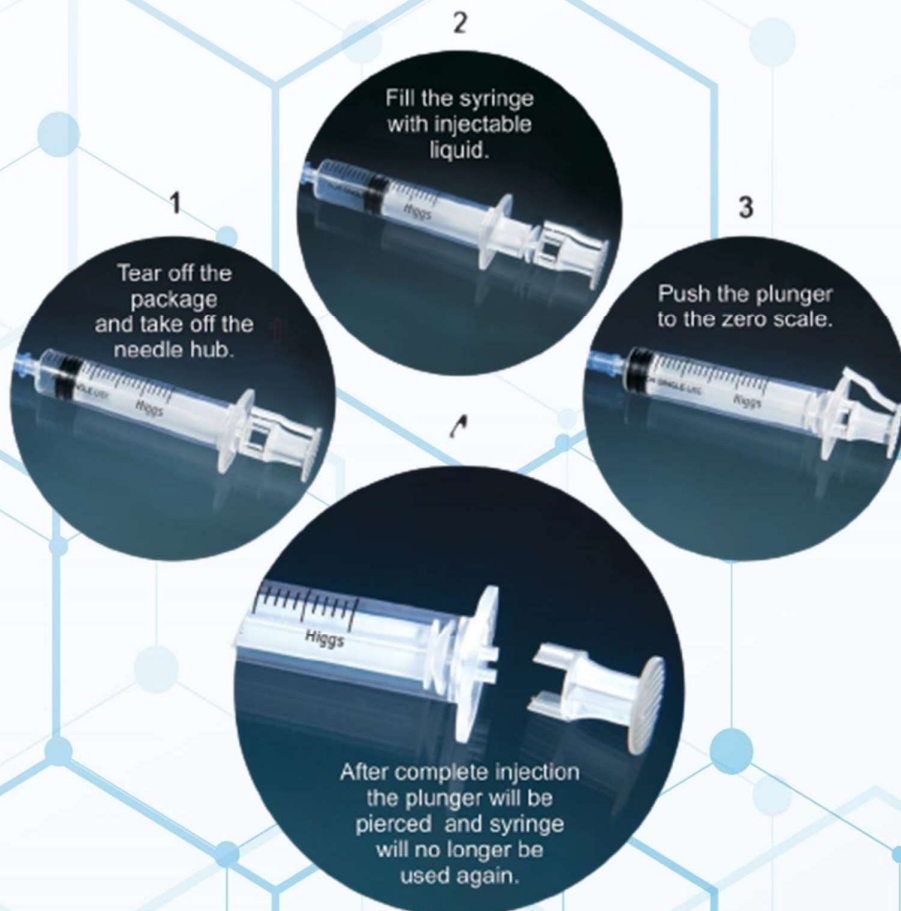

PRODUCT SPECIFICATIONS:

- Auto-disable (AD) syringes automatically become disabled after one use.
- AD Syringes are available with needle & without needle.
- Manufactured from medical grade virgin polypropylene.
- Prominent graduation ensures dosage accuracy and avoids wastage of valuable medicine.
- Latex free & thermoplastic Elastomer gasket is inert that provides minimum friction during movement and prevents leakage and back flow.
- Sharp triple facetted bevelled tip needle for minimum pain and extremely smooth penetration.
- Sterile, Disposable (Single use only), Non-Pyrogenic, Latex Free & Individually packed.
- Manufactured as per IS: 10258 / ISO: 7886-1.

PRODUCT PACKING:

- Individually packing: Ribbon Pack in tough laminated film, highly resistant to puncture and tearing.
- Individually packing: Blister Pack in medical grade paper and film making it safe and suitable for long-term storage.
- Inner packing: Corrugated carton box
- Outer packing: Corrugated carton box
- For more packaging details please download e-catalogue.

LOCK (DISABLE SYSTEM):

- Presence of loop & lock mechanism which automatically breaks the plunger, if reuse is attempted.

PRECAUTION:

- Before use, check the integrity of the sealing of the pack, if the pack is previously opened or damaged then do not use.
- Do not push the plunger forward before withdrawing medicament from the vial otherwise plunger will get locked.

WARNING:

- Do not attempt to re-sterilize and reuse.

PRODUCT RANGE:

Capacity	Types of syringes	Model of syringe	Standard Needle Dimensions
0.5 ml	3 Part Syringes	luer slip / luer lock	Needle: 24G, 26G, 28G, 30G, 31G
1 ml	3 Part Syringes	luer slip / luer lock	Needle: 24G, 26G, 28G, 30G, 31G
2 ml	3 Part Syringes	luer slip / luer lock	Needle: 18G, 20G, 22G, 23G, 24G
3 ml	3 Part Syringes	luer slip / luer lock	Needle: 22G, 23G, 24G
5 ml	3 Part Syringes	luer slip / luer lock	Needle: 18G, 20G, 21G, 22G, 23G, 24G
10 ml	3 Part Syringes	luer slip	Needle: 16G, 18G, 20G, 21G, 22G, 23G

Advantages:

- Note: As per customer requirement, we will provide any gauge and any length (0.5", 1", 1.25" & 1.5") of needles.



Our Prestigious Brand:



M/s. Hervik Healthcare Private Limited

Manufacturers of Disposable Medical Devices.

Reg. Office: Door Number. 8-202, Jamandla Street,
Chandragiri, Chittoor (Dist),
Andhra Pradesh - 517101, India.

Factory Address: Survey No: 600/B,
NH-140, Pichinaidupalli Road,
Chandragiri, Chittoor (Dist),
Andhra Pradesh - 517101, India.

Contact Mobile Number: +91 9985 979 770
Email: hervikhealthcare@gmail.com
Website: www.hervikhealthcare.com